



BladeCASSETTE - Single-Handed Scalpel Blade Remover

A simply safer way of removing blades in operating rooms.

Clicksmart BladeCASSETTE provides a simply safer way to remove and contain used scalpel blades in sterile operating environments. This world-first single-handed scalpel blade remover prevents the dangerous injuries commonly caused by removing blades with fingers, forceps or re-sheathing techniques.

- Wide base with grips to reduce slippage.
- Gamma radiation sterilisation.
- Transparent cartridges, easy counting of blades.
- Can be used in sterile packs.

With the single-use BladeCASSETTE you can simply 'push and click' to change and contain blades during an operation without risking the other hand. BladeCASSETTE stands and grips firmly on the drapes, making single-handed removal of blades possible. The other hand is out of harm's way until the blade is safely contained in the remover.

The BladeCASSETTE simple 'push and click' blade removal solution for steel scalpels:

- IMPROVES PATIENT SAFETY by enabling surgeons to use ergonomic precision steel scalpels.
- IMPROVES STAFF SAFETY by ensuring the used scalpel blade is removed with a safe single-handed technique and then safely contained until final disposal.


BladeCASSETTE can be included in sterile packs (e.g. custom kits or theatre packs), ensuring that the blade remover is immediately available for use when required. BladeCASSETTE can be sterilised by ETO or gamma radiation.

Counting of blades

BladeCASSETTE has been specifically designed to enable easy and accurate counting of blades after use in theatre.



Nodules help stabilise cassette by gripping drape.

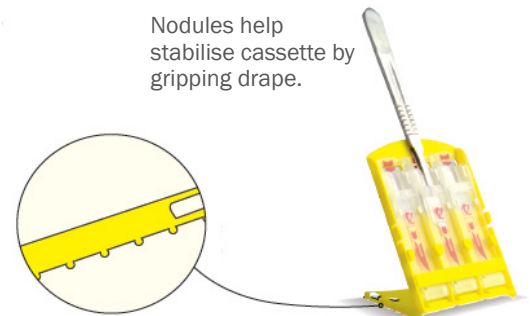

SINGLE-HANDED

Scalpel Blade Remover:
Safe & Easy to use



TWO-HANDED

Scalpel Blade Remover:
Unsafe & Difficult to use



Product Code	Description	QTY
QSSVCAS-3	Blade Cassette, 3 in 1 Blade Remover, Sterile	Box 100