



## Surgicric


Cricothyrotomy set to secure the airway in “cannot intubate - cannot ventilate” situations.

- Clear set composition. No confusion about different techniques
- Atraumatic specially designed dilator and tracheal tube
- Adjustable flange on the tracheal tube
- Aspiration of air via dilator
- New locking mechanism between dilator and tracheal tube
- Unique packing system

### Surgicric III

#### Seldinger technique

The set Surgicric III contains a unique introducer needle. This offers the advantage of an optimal skin puncture with temporary oxygenation via the 15 mm connector and the implementation of the Seldinger technique

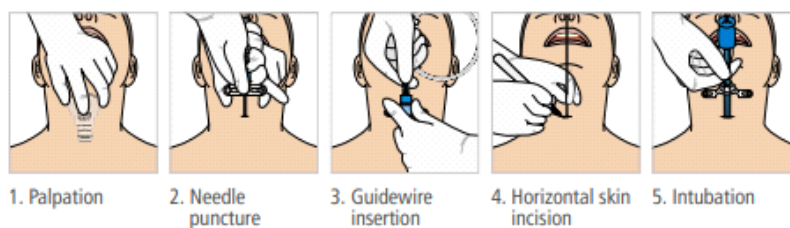

Description	QTY	Product Code
Surgicric III, includes safety scalpel, syringe 10ml, 6mm cuffed tracheal tube with dilator, extension tubing, necktape for tracheal tube, necktape for introducer needle, introducer needle & guidewire. Single use, sterile	Box 1	30-08-227-1