



Disposable Medical Proctoscopes

Diagnostic Line




The disposable proctoscope utilises a convenient battery-operated pen light, enabling examinations and minor therapeutic procedures to be performed in a non-specialist out-patients' department, on the ward or at home. This proctoscope can also be used with a plug-in light source produced by SapiMed, if desired. Proctoscopes are made of transparent plastic material and supplied in individually sealed bags. A free pen-light is included in each box of scopes.

TECHNICAL FEATURES


- Clear transparent plastic anosopes
- Range of four size/designs to suit varying clinical needs
- Flash free
- Guaranteed atraumatic
- Illumination by either clinic light source or pen-light
- Conforms fully with essential requirements of European Council Directive 93/42/EU
- Latex Free



A.4021

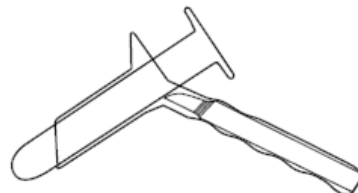

A.4018

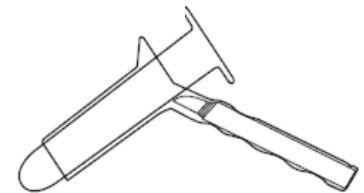

A.4019

Product Code	Description	QTY
A.4018	Proctoscope diagnostic/paediatric fluted beaked Disposable Tube length: 56mm, Tube diameter: 18mm For clear and easy rectal examination	Box of 75 pcs. + 1x pen-light
A.4019	Proctoscope diagnostic/paediatric trunk-end Disposable Tube length: 54mm, Tube diameter: 18mm For clear and easy rectal examination	Box of 75 pcs. + 1x pen-light
A.4021	"Trunk" operative anoscope, disposable. Tube length: mm 70 Tube diameter: mm 21	Box 25



Disposable Medical Proctoscopes

**Surgical
Line**




A.4022 - It has an increase of about 14mm over the A.4024, so offering the practitioner improved access higher up the rectum. The introducer is designed for easy and atraumatic introduction to the patient. A.4022 has similar design characteristics to code A.4023. It allows excellent control of the anal canal walls.

A.4023 - Flute beak to perform biopsies, elastic ligatures and haemorrhoid sclerosis, incision of intra-anal haematomas, excision of hypertrophic anal papillae, sutures with haemostatic stitches following polypectomy, I.R. coagulation etc.


A.4024 - A.4024 has similar design characteristics to code A.4023. Allows excellent control of the anal canal walls.



A.4022

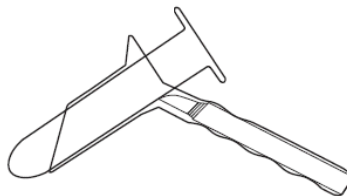

A.4023

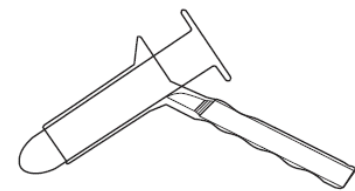

A.4024

Product Code	Description	QTY
A.4022	Surgical/Examination Blu Scope Proctoscope trunk-end Disposable Tube length: 79mm, Tube diameter: 23mm	Box of 50 pcs. + 1x pen-light
A.4023	Surgical/Examination Proctoscope flute beaked Disposable Tube length: 88mm maximum, Tube diameter: 23mm	Box of 50 pcs. + 1x pen-light
A.4024	Surgical/Examination Proctoscope trunk-end Disposable Tube length: 65mm, Tube diameter: 23mm	Box of 50 pcs. + 1x pen-light