




ScalpelCric

Cricothyrotomy set acc. to the scalpel technique

The ScalpelCric from VBM is a surgical set for cricothyrotomy according to the scalpel technique. The set has been developed to match the latest recommendations from the DAS 2015 guidelines*. Based on the simple description of 'stab, twist, bougie, tube', the ScalpelCric is meant to facilitate the execution of the scalpel technique.

Each component is designed to enhance the safety and efficacy of the technique. The set is packed in a specific sterile pouch with each component organised in the order of use. This reduces the risk of confusion at the critical time of the procedure. Scalpel #10, 14FR bougie (40 cm) with angled and rounded tip, 6.0 mm cuffed tube with smooth tip, extension tubing with suction capability, syringe for cuff inflation and necktape for tube fixation.

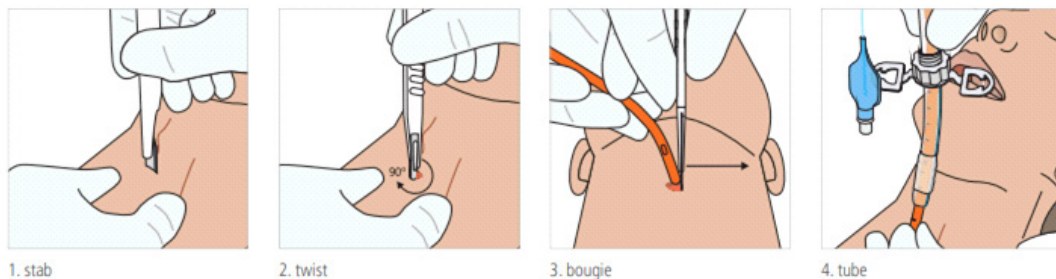

Set Includes

- Scalpel # 10
- 14FR Bougie (40 cm)
- 6.0 mm lockable tube
- Connecting hose
- Syringes 10 ml
- Collar

Description	QTY	Product Code
ScalpelCric, sterile	Box 1	30-08-447-1
ScalpelCric, non sterile (for training)	Box 1	30-08-447